


BTK (Ab-223) Antibody


Catalog No.	75675
Product Name	BTK (Ab-223) Antibody
Size	50ug, 100ug
Concentration	1 mg/ml
Specificity	Detects endogenous levels of total BTK protein.
Reactivity	Human, Mouse
Immunogen	Synthesized non-phosphopeptide derived from human BTK around the phosphorylation site of tyrosine 223 (A-L-Y ^P -D-Y).
Applications	IHC: 1/50 – 1/100; WB: 1/500 – 1/1000
Swiss-Prot No.	Q06187 (Human)
Host Species	Rabbit
Clonality	Polyclonal
Isotype	IgG
Storage	-20°C for one year. Avoid freeze / thaw cycles.
Purification	Affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.

Images

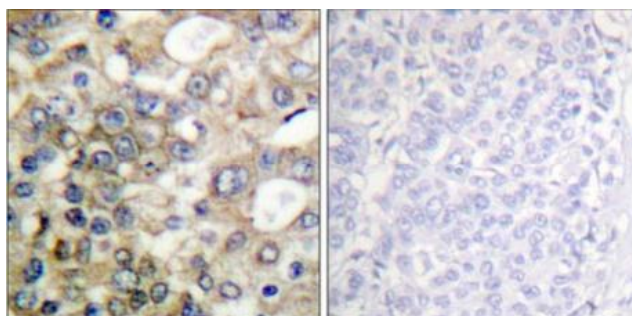

Immunohistochemical analysis of paraffin-embedded human breast carcinoma tissue using BTK (Ab-223) antibody.

Left: peptide –
Right: peptide +

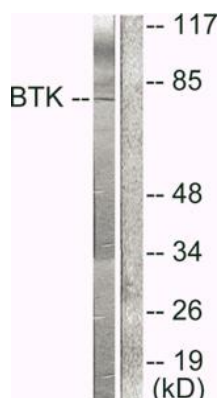

Western blot analysis of extracts from Hela cells, using BTK (Ab-223) antibody.

Left lane: peptide –
Right lane: peptide +

For Research Use Only